

Quick User Guide to Comfort Audio FM system - Body worn TOD version

Components

Teachers
Microphone



Pupils receiver with
neckloop



STEP 1: Switch on teachers microphone


Press and hold the top button and look at the screen on the back of the unit. Press until the GREEN bar has fully illuminated on the screen.



A STEP 2: Fitting and switching on pupils receiver

Place the neckloop around hearing aid users neck and plug in to receiver as shown in diagram. Switch hearing aid to the T Position.


As with the microphone press and hold the top button to switch on.

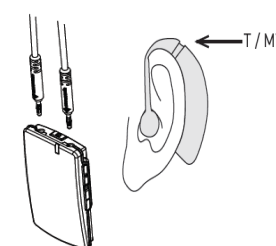

STEP 3: Wearing the microphone and receiver

Clip on microphone to clothing, the microphone works best around 30cm.

A necklace is provided if the clip is not possible to use.

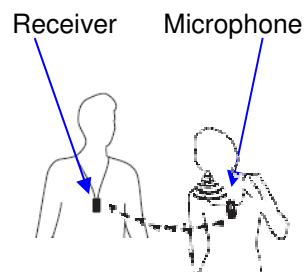

STEP 4: Increasing the volume


The volume of the neckloop can altered via the + and - switch on the side of the receiver.



STEP 5: Charging


Plug charger in to the mains using the charger provided. The indicator will be RED when charging.

Once fully charged the indicator will go GREEN.

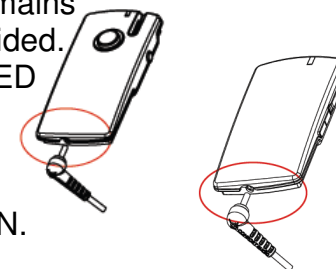

Connecting an audio source


An audio source can be connected e.g. White board audio. Simply find the OUTPUT from the device. This will say LINE OUT or HEADPHONE SOCKET. Use connection lead provided and connect to Audio input of the microphone.

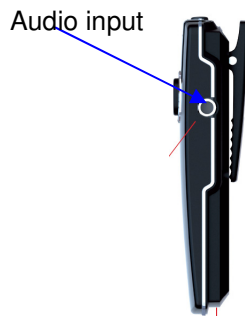

Changing the channel (if required)

If there is interference the channel may require changing. Push jog wheel (shown below)

In menu select *communication* select *channel*. Select channel with a green dot.

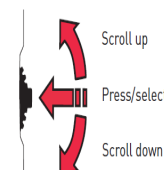

Pairing the system (if required)

After the channel has been changed and there is still interference the system may require pairing.

Press and hold the receivers button for 2 seconds. Double click the pairing button on the transmitter.

